BVS
BEYOND VALUABLE SOLUTION

cobetter®
filtration

Biophsep™
Hollow Fiber



www.bvskorea.com
02-548-6211



Product Features

- High flow rate, high filtration capacity
- Modified hydrophilic PES hollow fiber membrane provides low protein adsorption, less membrane fouling and easy cleaning
- Integrated device without additional assembly or device holder, quick installation and operation
- Re-washable with 0.5M NaOH solution
- Simple and reliable linear scale-up

Cobetter Biophsep™ hollow fiber with good filtration resolution and low fouling , helps biopharmaceutical customers to improve product yield and filtration load, and reduce production costs.

Typical Use by Application

- Purification, concentration and diafiltration of vaccine
- Purification, concentration and diafiltration of viral vector
- Clarification of cells and bacterial in fermentation broth
- Recycling and washing of cells and bacteria
- Concentration and diafiltration of protein



Hollow Fiber Filter (Autoclavable) Part Number Description(1.0mm)

Module	Effective Length (cm)	Process Volume	MWCO	Fiber Inner Diameter (mm)	Effective Area (cm^2/m^2)	Fiber Number	Flow rate @2000s^{-1} (mL/min)	Flow rate @6000s^{-1} (mL/min)	Overall Dimension (mm*mm)	Inlet/ Outlet	Side Port
Mini	30	<250ml	750R=750kD-R M20R=0.2um-R	1.0	28	3	35	106	10*346	Luer Connection Female	Luer Connection Female
	60	<500ml		1.0	56	3	35	106	10*646		
Minilab	30	<750mL		1.0	75	8	94	282	10*346		
	60	<1.5L		1.0	150	8	94	282	10*646		
Lab	30	<1.5L		1.0	170	18	212	635	13*348	0.5'' TC (Disk Outer Diameter 25mm)	3/16'' HB
	60	<3L		1.0	340	18	212	635	13*648		
Lab+	30	<5L		1.0	0.05	60	706	2118	18*355	0.5'' TC (Disk Outer Diameter 25mm)	3/8'' HB
	60	<10L		1.0	0.10	60	706	2118	18*655		
Pilot	30	<10L		1.0	0.12	130	1529	4588	33*361	1.5'' TC (Disk Outer Diameter 50.5mm)	0.5'' TC (Disk Outer Diameter 25mm)
	60	<20L		1.0	0.25	130	1529	4588	33*661		
	110	<40L		1.0	0.45	130	1529	4588	33*1161		
Pilot+	30	<15L		1.0	0.15	160	1882	5647	33*361		
	60	<30L		1.0	0.30	160	1882	5647	33*661		
	110	<50L		1.0	0.55	160	1882	5647	33*1161		
Column Max	60	<400L		1.0	4.2	2550	30000	90000	101.5*622	2.0''TC (Disk Outer Diameter 64mm) For use with SS housing	1.5''TC (Disk Outer Diameter 50.5mm) For use with SS housing

HF-E-XX-XXXX-XX-XX-A-P e.g.:HFEMIM20R1060AP

Product type
HF Hollow Fiber

Membrane
E PES

Module
MN Mini
MI Minilab
LA Lab
LP Lab+
PI Pilot
PP Pilot+
CM Columm Max

MWCO
750R 750kD-R
M20R 0.20μm-R

Fiber Inner Lumen
10 1.0mm

Effective Length
30 30cm
60 60cm
11 110cm





Sterilization
A Autoclavable

P Pharmaceutical


Material Construction

Module Component	Material	Advantages and Functions
Hollow Fiber Membrane	mPES	The modified hydrophilic PES hollow fiber membrane features low adsorption, less membrane fouling and continuous high flow rate for faster processing time, the membrane is proved to effectively retain virus particles so as to achieve the purpose of concentration and buffer exchange.
Potting Glue	Polyurethane/Epoxy	Potting glue wraps each hollow fiber to provide support for the hollow fiber membrane, and at the same time completely separates the inlet flow channel and the permeate flow channel.
End Cap	White Polysulfone	Flow channel connection for liquid in and out, with good chemical compatibility.
Shell	Transparent Polysulfone	Connect the inlet and outlet caps to form an integrated assembly, while providing a flow channel for the permeate, with good chemical compatibility.


Overall Dimension

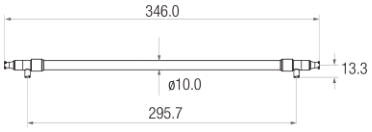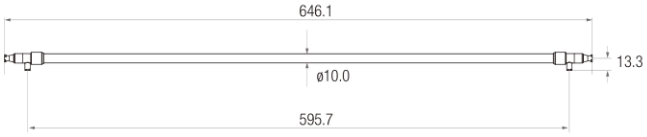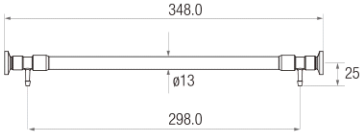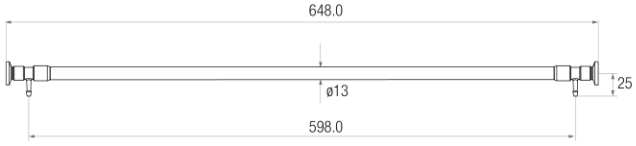	MI-Minilab -30cm
	MI-Minilab-60cm
	Lab-30cm
	Lab-60cm


Overall Dimension

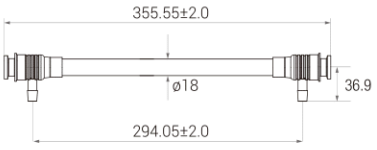


Lab⁺-30cm

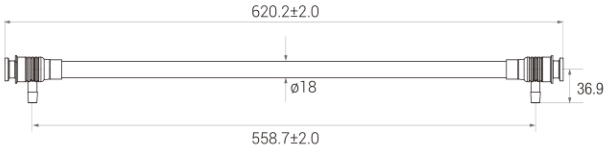


Lab⁺-60cm

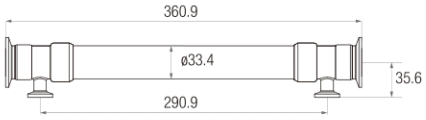

Pilot / Pilot⁺ 30cm

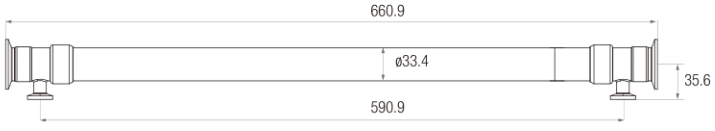


Pilot / Pilot⁺ 60cm

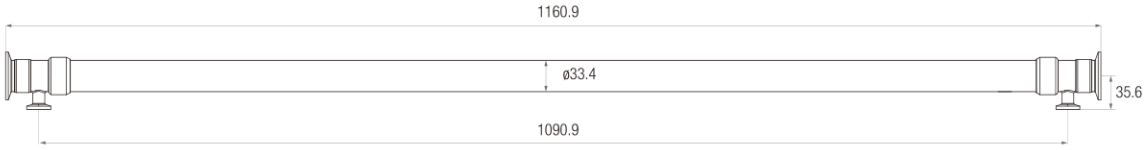


Pilot / Pilot⁺ 110cm


Overall Dimension

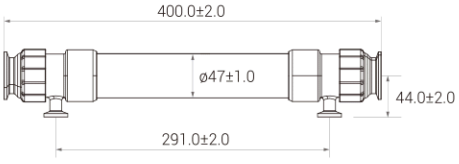


Mini-Process
30cm

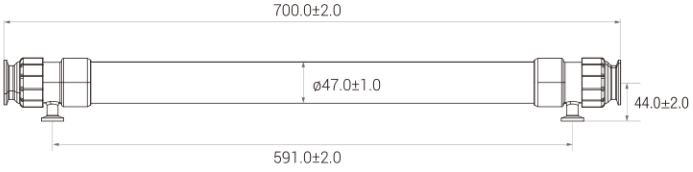


Mini-Process
60cm

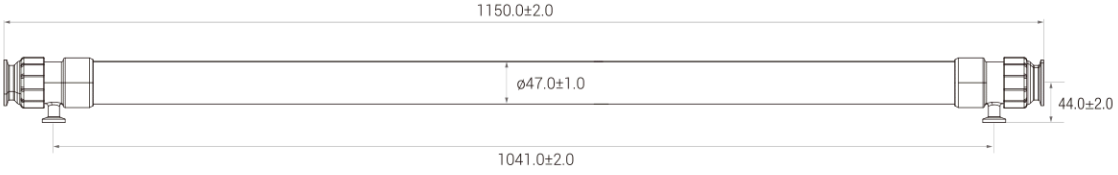


Mini-Process
110cm

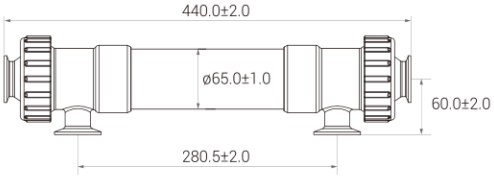

Mid-Process
30cm

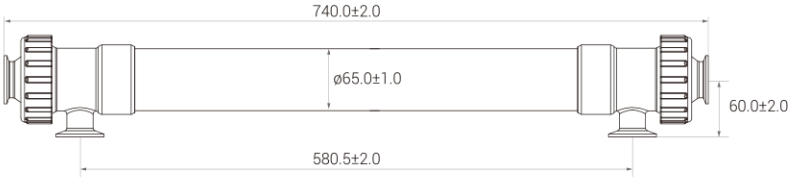


Mid-Process
60cm

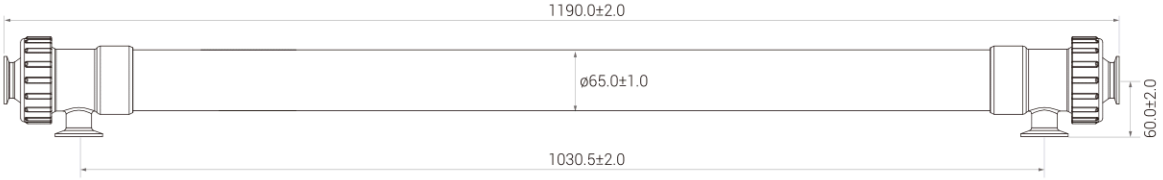


Mid-Process
110cm


Overall Dimension

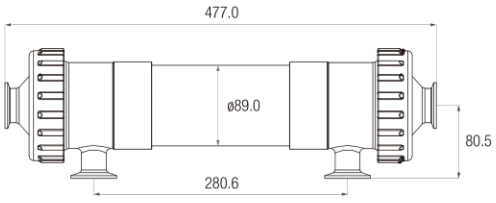


Process 30cm

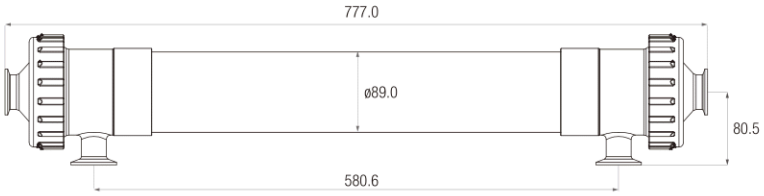


Process 60cm

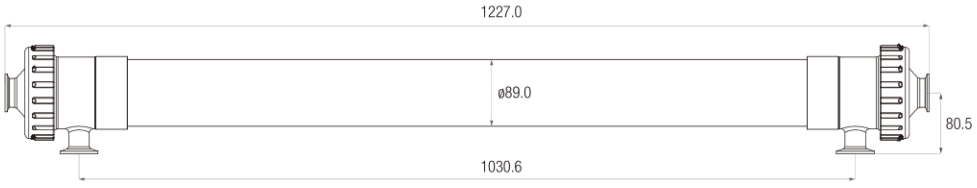

Process 110cm

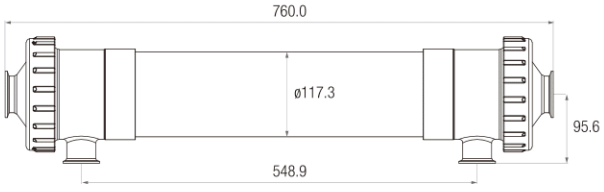

Maxi 60cm

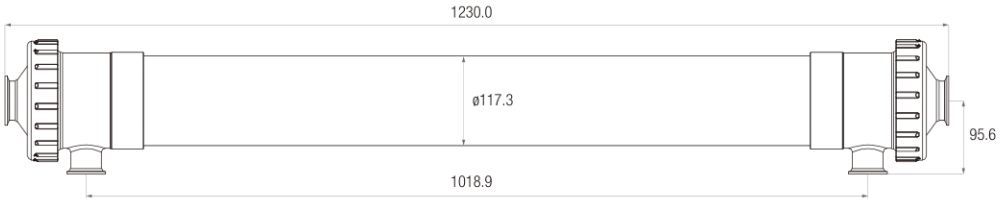

Max 110cm

Key Figures

Shear Rate

The circulating flow rate of the hollow fiber depends on the product tolerance to the shear rate, in generally it will be set at 2000/s to 10000/s, which is much smaller than the shear rate generated by the turbulent flow on the surface of the cassette screen.

For general materials, we usually choose a shear rate of 4000/s to 6000/s;

If the product is sensitive to shear force (such as lentivirus, coronavirus or macromolecular protein expressed by animal cells, large plasmids and LNP, etc.), the shear rate needs to be reduced to 2000/s;

If the product is with good resistance to shear force (such as small molecular proteins expressed by bacteria, etc.) shear rate could be increased to 8000/s to 10000/s.

Fiber Inner Lumen

Module with fiber of 1.0mm ID is ideal for product with high cell density or high solid content or high viscosity.

Module with fiber of 0.5mm ID is widely used in most application scenarios to improve mass transfer efficiency.

Only module with fiber of 1.0mm ID is autoclavable, as steam cannot enter the fiber of 0.5mm ID effectively.

MWCO

Sufficient separation efficiency and the risks of fouling should both be taken into consideration when choosing MWCO.

To minimize the risk of fouling that inevitably occurs during membrane processing, it is advisable to choose membranes with relatively smaller pore sizes, as long as the separation efficiency and flux are maintained. This helps reduce the intrusion of impurity particles into the membrane pores, thus extending the membrane's service life.

Common processing scenarios are as following:
· Heat source removal of small molecule in Chinese medicines : 5kD, 10kD
· Ultrafiltration and concentration of mABs, buffer exchange : 30kD, 50kD
· Concentration, purification and removal of virus : 100kD, 300kD, 500kD, 750kD
· Clarification of recombinant protein, antibody : 500kD, 750kD
· Concentration for bacteria: 500kD, 750kD
· Clarification of lysate: 0.2μm, 0.45μm, 0.65μm

Effective Length

The process scale-up feature of hollow fibers is that: direct process scale-up can be carried out as long as the effective length is kept the same.

On the opposite, due to the significant pressure drop difference between inlet and outlet existed in different lengths, the internal pressure and flow velocity distribution of the flow channel also change correspondingly, so linear scale-up can not be performed on devices with different length. When processing with more fouling and high viscosity product, it is preferred to choose module with shorter length.



www.bvskorea.com
02-546-6211
서울시 서초구 나루터로 70 영서빌딩 엠피스센터 111호